RadiForce®

GX530-BL GX530-CL GX530-CLAR

Monochrome LCD Monitor

English

Deutsch

Français

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Instructions for Use: ■ Setup Manual

Important: Please read PRECAUTIONS, this Setup Manual and the User's Manual stored on the CD-ROM carefully to familiarize yourself with safe and effective usage. Please retain this manual for future reference.

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Compatible Resolutions

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The monitor supports the following resolutions.

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^{*1} A supported graphics board is required to display.

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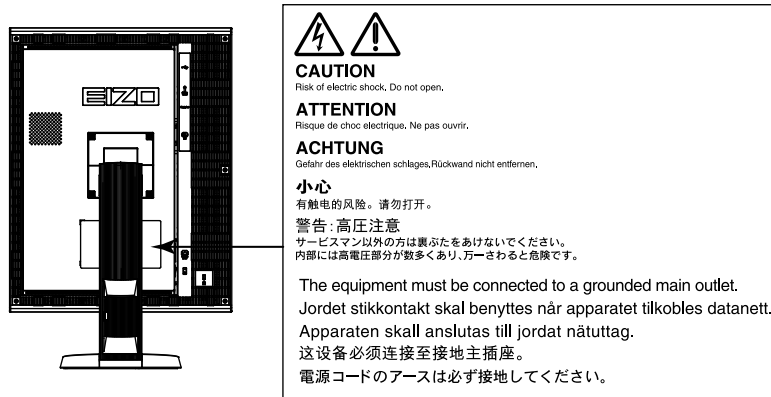
RadiForce® GX530-BL GX530-CL GX530-CLAR

Monochrome LCD Monitor

Location of the Caution Statements

Symbols on the unit

Symbol	This symbol indicates
	Main Power Switch: Press to turn the monitor's main power off.
	Main Power Switch: Press to turn the monitor's main power on.
	Power button: Press to turn the monitor's power on or off.
	Alternating current
	Alerting electrical hazard
	CAUTION: Refer to SAFETY SYMBOLS section in PRECAUTIONS.

No part of this manual may be reproduced, stored in a retrieval system, or transmitted, in any form or by any means, electronic, mechanical, or otherwise, without the prior written permission of EIZO Corporation. EIZO Corporation is under no obligation to hold any submitted material or information confidential unless prior arrangements are made pursuant to EIZO Corporation's receipt of said information. Although every effort has been made to ensure that this manual provides up-to-date information, please note that EIZO monitor specifications are subject to change without notice.


Package Contents

Check that all the following items are included in the packaging box. If any items are missing or damaged, contact your local dealer.


Note

- Please keep the packaging box and materials for future movement or transport of the monitor.


- Monitor
- Power cord

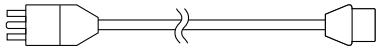
- Digital signal cable (DisplayPort) : PP300

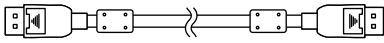- Digital signal cable (Dual Link) : DD300DL

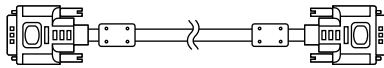
- USB cable: UU300

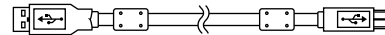- EIZO LCD Utility Disk (CD-ROM)
- Setup Manual (this manual)
- PRECAUTIONS
- Cleaning Kit "ScreenCleaner" (for GX530-CLAR only)

Controls and Functions

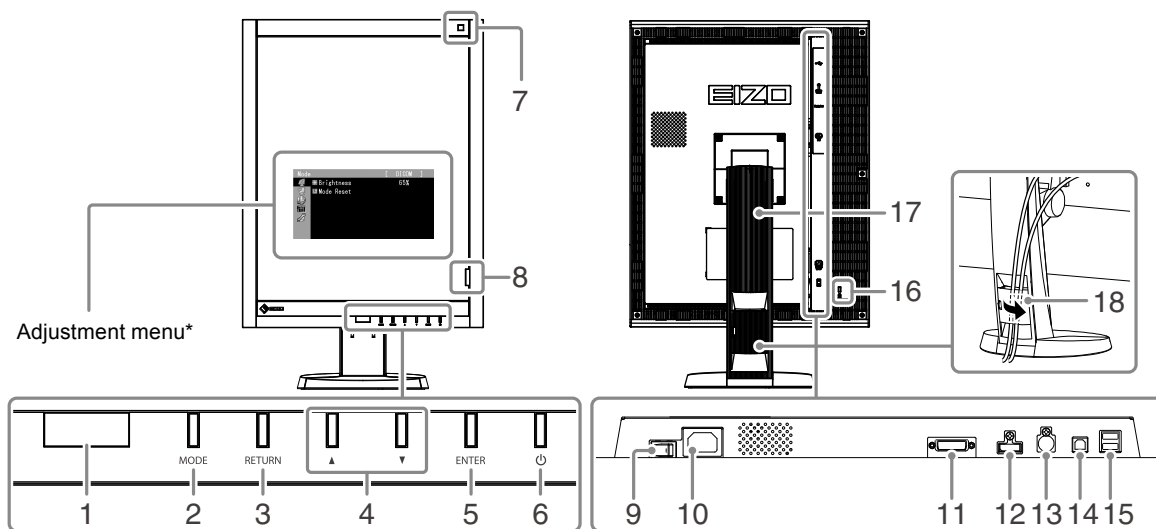1. Presence Sensor	Detects the movement of a person in front of the monitor.
2. MODE button	Switches the CAL Switch mode.
3. RETURN button	Cancels the setting/adjustment and exits the Adjustment menu.
4. ▲▼ button	Chooses and increases/decreases adjustment settings for advanced adjustment using the Adjustment menu.
5. ENTER button	Displays the Adjustment menu, determines an item on the menu screen, and saves values adjusted.
6. ⏻ button	<ul style="list-style-type: none"> • Turns the power on or off. • Indicates the monitor's operation status. Green: Operating, Orange: Power saving mode, Off: Main power/power off
7. Illuminance sensor	Measures the environmental illumination.
8. Integrated Front Sensor	Used for calibration and Grayscale Check.
9. Main power switch	Turns the main power on or off.
10. Power connector	Connects the power cord.
11. Input signal connector	DVI-D connector
12. Input signal connector	DisplayPort connector
13. PS/2 connector	Used for maintenance. Do not use for other applications. This may cause the monitor to malfunction.
14. USB port (UP)	Connects the USB cable to use the software that needs USB connection, or to use USB Hub function.
15. USB port (DOWN)	Connects a peripheral USB device.
16. Security lock slot	Complies with Kensington's MicroSaver security system.
17. Stand	Used to adjust the height and angle of the monitor screen.
18. Cable holder	Covers the monitor cables.

* For instructions for use, refer to the User's Manual (on the CD-ROM).

Connecting Cables

Attention

- Check that the monitor and the PC are powered off.
- When replacing the current monitor with this monitor, be sure to change the PC settings for resolution and vertical scan frequency to those that are available for this monitor. Refer to the compatible resolution table (back of cover page) before connecting the PC.

Note

- When connecting multiple PCs to this product, refer to the User's Manual (on the CD-ROM).


1. Turn the monitor screen 90° in clockwise direction.

The landscape position is the default monitor orientation. Turn the monitor screen ninety degrees to portrait position before installing it.

2. Connect the signal cables to the input signal connectors and PC.


Check the shape of the connectors, and connect the cables.

After connecting the signal cable, tighten the screws of the connectors to secure the coupling.

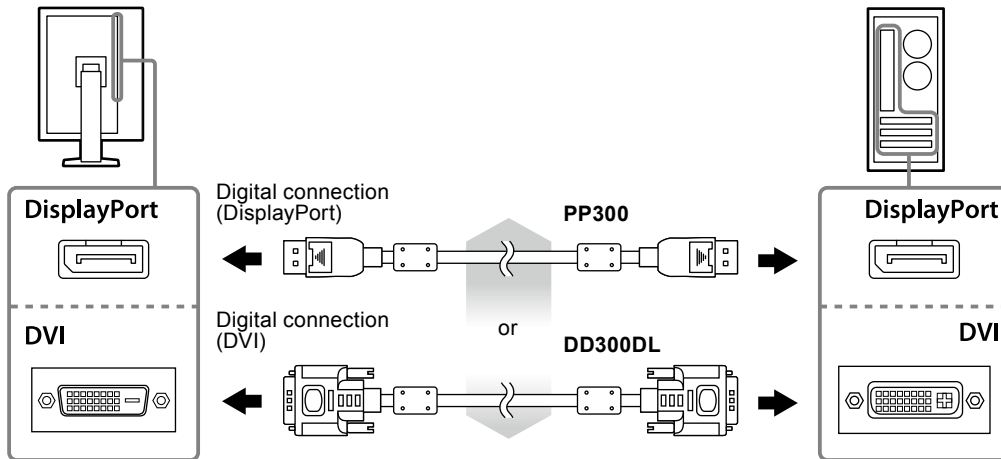3. Plug the power cord into a power outlet and the power connector on the monitor.

4. Connect the USB cable when using RadiCS LE or ScreenManager Pro for Medical.

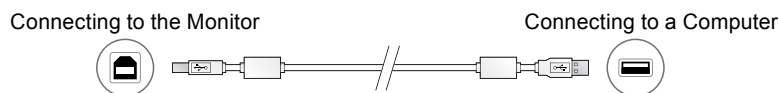5. Press \mathcal{P} to turn on the monitor.

The monitor's power indicator lights up green.

6. Turn on the PC.

The screen image appears.


If an image does not appear, refer to "No picture" (page 6) for additional advice.

Attention

- Turn off the monitor and PC after using them.
- For the maximum power saving, it is recommended that the Power button be turned off. Turning off the main power switch or unplugging the power cord completely shuts off power supply to the monitor.

Adjusting the Screen Height and Angle

Hold left and right edge of the monitor with both hands, and adjust the screen height, tilt and swivel of the screen to the best condition for working.

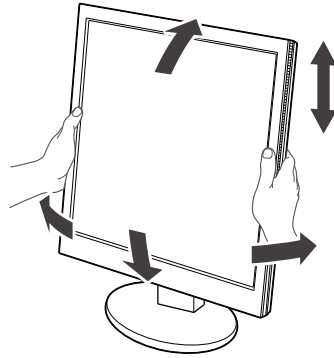
Attention

- Be sure that the cables are correctly connected.
-

Showing Button Guide

Press the front buttons (except  button), the Button Guide appears above the button.

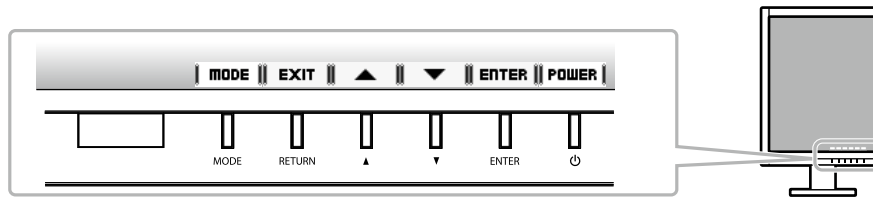
Note

- The Button Guide will continue to appear while the Adjustment menu or mode menu is showing.
 - The Button Guide is displayed differently depending on the selected menu or status.
-

Selecting Display Mode


This function allows easy selection of a desired mode according to monitor application.

Mode	Purpose
1-DICOM	Available to display in DICOM mode.
2-Custom	Available for the setting according to your preference.
3-CAL	Displays the screen adjusted by calibration software.

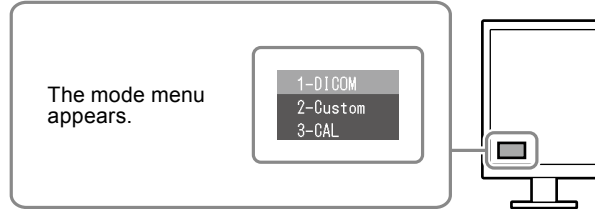Procedure

1. Press MODE.
The mode menu appears at the bottom left of the screen.
2. The mode among the list is highlighted in turn each time MODE is pressed.
You can switch the mode with ▲ or ▼ while the mode menu is displayed.

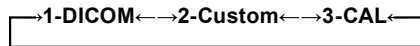Controlling Monitor Quality

The quality control software “RadiCS LE” (for Windows) is included with this monitor.

To calibrate the monitor and manage the calibration history, install RadiCS LE with reference to the RadiCS LE User's Manual (on the CD-ROM). Use the RadiCS SelfQC function and built-in Integrated Front Sensor to perform periodic calibration and Grayscale Check for the monitor independently. For details, refer to the User's Manual (on the CD-ROM).

Attention

- Connect the monitor to the PC with a USB cable when using RadiCS LE.

Note

- To perform advanced quality control that conforms to medical standards, use the optional monitor quality control software kit “RadiCS UX1”.

No picture


If a problem still remains after applying the suggested remedies, contact your local dealer.

Problem	Possible cause and remedy
<p>1. No picture</p> <ul style="list-style-type: none"> Power indicator does not light. <hr/> <ul style="list-style-type: none"> Power indicator is lighting orange or green. <hr/> <ul style="list-style-type: none"> Power indicator flashes orange and green. 	<ul style="list-style-type: none"> Check whether the power cord is connected properly. Turn the main power switch on. Press . Turn off the main power, and then turn it on again a few minutes later. <hr/> <ul style="list-style-type: none"> Switch the input signal. For details, refer to the User's Manual (on the CD-ROM). Operate the mouse or keyboard. Check whether the PC is turned on. When the Presence Sensor is set at a value from "Level-1" to "Level-4", the monitor may have shifted to the power saving mode. Come near to the monitor. <hr/> <ul style="list-style-type: none"> There is a problem in the device that uses DisplayPort connection. Solve the problem, turn off the monitor, and then turn it on again. Refer to the user's manual of the device that uses DisplayPort connection for details.
<p>2. The message below appears.</p> <ul style="list-style-type: none"> This message appears when no signal is input. Example: <div data-bbox="316 949 627 1149" data-label="Image"> <p>A screenshot of a black error message box. At the top, it says "No Signal" in white. Below that is a small icon of a DisplayPort connector and the text "DisplayPort". Underneath, it lists "fH: 0.0kHz" and "fV: 0.0Hz". At the bottom, it says "Check the input signal." in white.</p> </div> <ul style="list-style-type: none"> The message shows that the input signal is out of the specified frequency range. Example: <div data-bbox="316 1261 627 1503" data-label="Image"> <p>A screenshot of a black error message box. At the top, it says "Signal Error" in white. Below that is a small icon of a DVI connector and the text "DVI Digital". Underneath, it lists "fD: 135.0MHz", "fH: 79.9kHz", and "fV: 75.0Hz". At the bottom, it says "Check the input signal." in white.</p> </div> <p>fD: Dot Clock fH: Horizontal Scan Frequency fV: Vertical scan frequency</p>	<p>This message appears when the signal is not input correctly even when the monitor functions properly.</p> <ul style="list-style-type: none"> The message shown left may appear, because some PCs do not output the signal soon after power-on. Check whether the PC is turned on. Check whether the signal cable is connected properly. Switch the input signal. For details, refer to the User's Manual (on the CD-ROM). <hr/> <ul style="list-style-type: none"> Check whether the PC is configured to meet the resolution and vertical scan frequency requirements of the monitor (see "Compatible Resolutions"). Reboot the PC. Select an appropriate display mode using the graphics board's utility. Refer to the manual of the graphics board for details.

RadiForce® GX530-BL GX530-CL GX530-CLAR

LCD-Monochrommonitor

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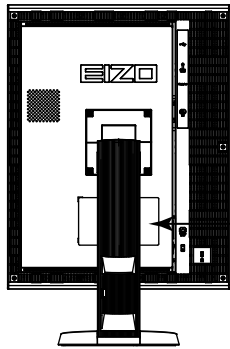CAUTION
Risk of electric shock. Do not open.

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